

# CIRCULAR STAPLER CURVED

## ORDER INFORMATION

ARTICLE NO.	FIGURE	DESCRIPTION	COLOUR CODE	STAPLE QUANTITY	ANVIL OUTSIDE DIAMETER	LUMEN SIZE
CSC21		Conventional	Orange	16	21.0 mm	13.0 mm
CSC21-A		Air-Sealed	Orange	16	21.0 mm	13.0 mm
CSC25		Conventional	White	20	25.0 mm	17.0 mm
CSC25-A		Air-Sealed	White	20	25.0 mm	17.0 mm
CSC29		Conventional	Blue	24	29.0 mm	20.5 mm
CSC29-A		Air-Sealed	Blue	24	29.0 mm	20.5 mm
CSC33		Conventional	Green	28	33.0 mm	24.8 mm
CSC33-A		Air-Sealed	Green	28	33.0 mm	24.8 mm

Manufacturer acc. to directive 93/42/EEC:  
 Touchstone International Medical Science Co., Ltd. | 21 A Science Plaza  
 International Science Park | No. 1355 Jinjihu Avenue | Suzhou | China

Distributed by:  
 B. Braun Surgical, S.A. | Carretera de Terrassa, 121 | 08191 Rubí | Spain  
 Phone +34 93 586 6200 | Fax +34 93 588 1096 | [www.bbraun.com](http://www.bbraun.com)

## AESCULAP® – a B. Braun brand

Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany  
 Phone +49 7461 95-0 | Fax +49 7461 95-2600 | [www.aesculap.com](http://www.aesculap.com)

The main product trademark "Aesculap" is a registered trademark of Aesculap AG.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.



GENERAL  
 OPEN SURGERY

## CSC – CIRCULAR STAPLER CURVED




SINGLE USE INTERNAL STAPLER

# CIRCULAR STAPLER CURVED

SINGLE USE INTERNAL STAPLER


## PRODUCT FEATURES

- RBS – Rotation Blade System. Easy & precise cutting.
- ASL – Built-in Automatic Safety Lock. Prevents accidental mis-firing.
- Bendos – Feeding tube re-entry pulley, protects anastomotic orifice.
- PSS – Petal Shaped Shaft. Easy insertion with firm attachment.
- Adjustable staple closure height ranging from 1 mm to 2.5 mm.

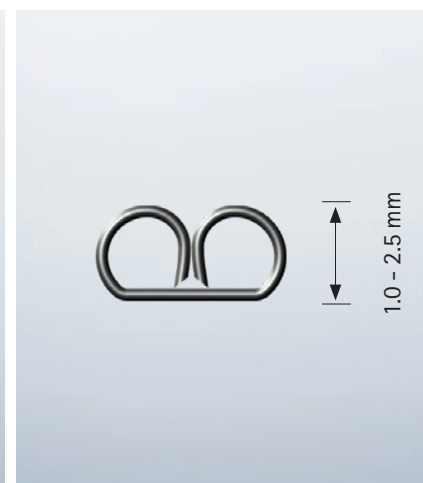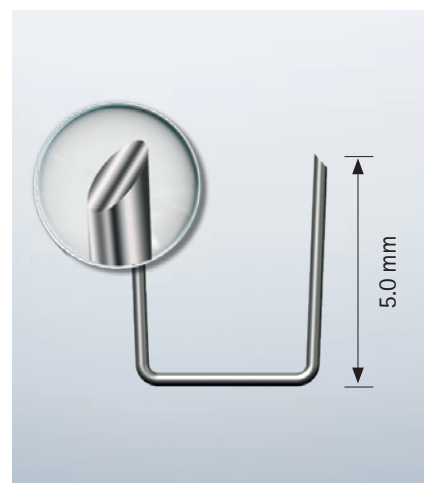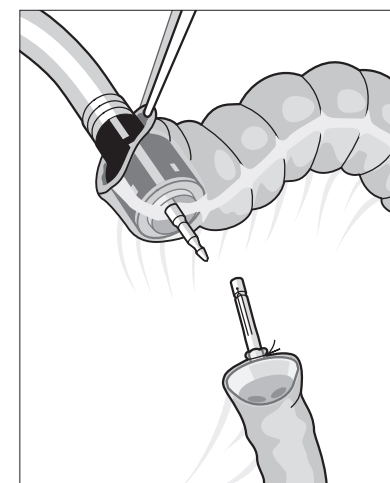

Before firing

After firing

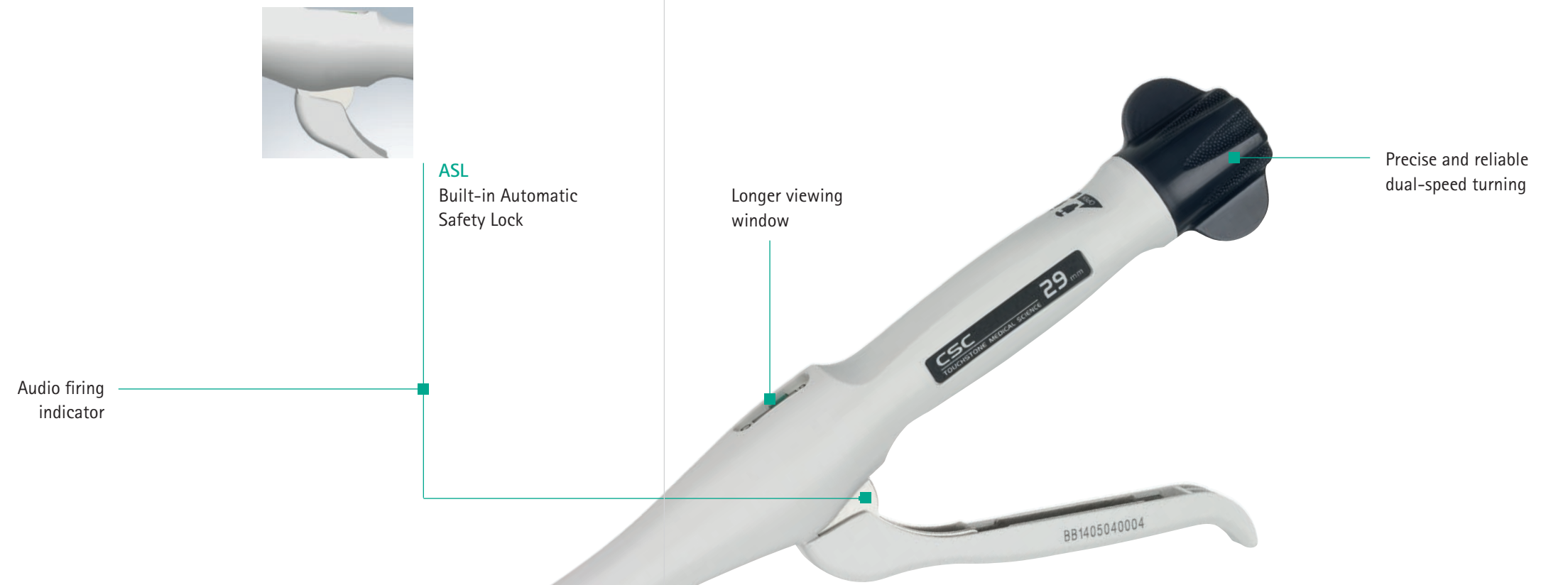

Audio firing indicator

ASL  
Built-in Automatic Safety Lock


Longer viewing window

Precise and reliable dual-speed turning


## DESCRIPTION



Bendos  
Feeding tube re-entry pulley

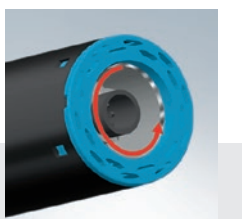

RBS Easy-cut  
Rotation Blade System

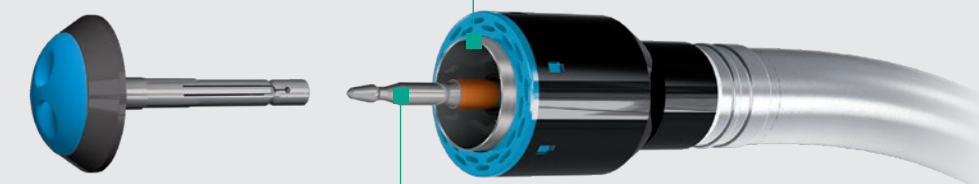

Slim trocar for efficient piercing